

INVACARE®

# RELIANT™ 600 HEAVY DUTY POWER LIFT

WITH MANUAL LOW BASE



Designed to meet the special needs of your residents weighing up to 600 lbs., the **Invacare® Reliant™ 600 Heavy Duty Power Lift with Manual Low Base** is secure and sturdy, providing a safe alternative to manual lifting. Both durable and functional, the Reliant 600 Lift is also a reassuring and comfortable resident lift. When answering the needs of every resident is your priority, this is your wisest choice.



**Yes, you can.®**

---

## **Invacare® Reliant™ 600 Heavy Duty Power Lift with Manual Low Base Features**

High lift range allows transfer to high surfaces

Low friction casters significantly improve rollability and contribute to caregiver safety

Designed to meet the needs of the growing bariatric population up to 600 lb.

Manual lowering for instant response to power loss

Unique sling fabric designed to provide the utmost in comfort and strength

Available with power opening base (model # RPL600-2)

Essential device for safe resident handling in hospitals and LTC (long-term care) facilities

Limited warranty includes three years on lift, one year on electric components

Rear locking casters are standard

Manufactured in an ISO13485 facility

---

Product Code: RPL600-1




RPL600-1: Powered with low base



RPL600-2: Powered with power opening low base

---

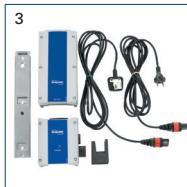

#### 1. **Manual Lowering**

Manual emergency lowering gently and safely lowers the resident in case of total power loss

Essential safety feature

Manual lowering is easy and effortless. Simply insert the end of a ball point pen into the hole marked "emergency" or lift on the accessible ring-pull

#### 2. **Invacare® Reliant™ Digital Scale**

Retrofits on all Reliant 450/600 lifts

Easy to operate, four-button display

Automatic shut-off feature prolongs battery life

Two-year limited warranty

#### 3. **Electronics**

Removable battery box conveniently located on mast features an emergency off button and beeps when battery needs recharged. Requires one 24-volt battery



Quick-release battery eliminates downtime

Anti-entrapment feature stops motion if boom meets resistance while lowering

#### 4. **Charging Station**





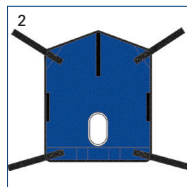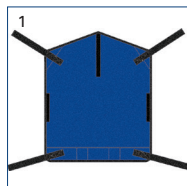







Battery can be replaced immediately with second battery in wall-mounted charger

UL listed and CSA certified



#### 1. **Mesh sling model # R140**

#### 2. **Mesh sling with commode opening model # R141**

							
Model #	Height At Sling Hook Up	Base Dimensions	Casters	Weight	Weight Capacity	Battery Life Indicator	Emergency Controls
RPL600-1	28"-68" H	4.5"H Clearance 41"W Open 26.5"W Closed 48"L	3.5" Front 5" Rear	112 lb.	600 lb.	N/A	Stop and Down
RPL600-2	28"-68" H	4.5"H Clearance 41"W Open 26.5"W Closed 48"L	3.5" Front 5" Rear	115 lb.	600 lb.	N/A	Stop and Down
							
	Battery	Charge Time	Charger Input	Charger Output	Limited Warranty		
RPL600-1	2-24V DC Rechargeable	6 Hours	100-240V AC	29.5 DC	3 Years Lift 1 Years Electronics		
RPL600-2	2-24V DC Rechargeable	6 Hours	100-240V AC	29.5 DC	3 Years Lift 1 Years Electronics		

**Invacare Thailand Ltd**

No. 399 Interchange 21,  
33rd Floor, Sukhumvit Road,  
North Klongtoey Sub-district,  
Wattana District,  
Bangkok Metropolis  
Phone: +66 0 2108 8555  
Web: [www.invacare.com](http://www.invacare.com)

©2013

Invacare Corporation.  
All rights reserved. Trademarks  
are identified by the symbols  
™ and ®. All trademarks are  
owned by or licensed to Invacare  
Corporation unless otherwise noted.  
Specifications are subject to change  
without notification.

Form No. 13-176  
130304

